



Portable  
Room-to-Room  
Unique Design



Our Roomer S overhead lift is a portable lift that offers the possibility to transfer a patient from one room to another in one single motion without making modifications to doorway headers.

#### A Comfortable Experience

The Roomer S offers optimum comfort when used with our sling portfolio. The room-to-room function reduces the number of lifts required and minimizes strain on both the patient and the caregiver. It offers a soft-start, soft-stop and constant lifting speed regardless of the patient weight and direction, making lifting a comfortable experience. The lift has a weight capacity of 220kg/484lb and can be used with either a 2, 3 or 4-point suspension.

#### Unique Design

The Roomer S is suitable for rooms with sliding doors or fire doors as no modification to the doorway is required, thus making installation

much easier. For facilities, its unique design eliminates the need for turntables and rail switches. The Roomer S also features a lightweight quick release hanger bar in two different widths to accommodate most lifting situations. Charging can be done through the hand control for optimal efficiency.

#### Lift Intelligence

The Roomer S features an advanced technology platform that allows tracking of lift usage and visualization of collected data. This integrated lift intelligence allows data to be retrieved using our Data and Service Application on a tablet or computer. Detailed information can be viewed directly in the application or exported to other media for additional analysis.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone: +46 8 665 35 00  
Fax: +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

# Accessories



**Hanger bar**  
Quick release for optimal convenience. Order separately. Part no.: 55751 (45cm/17.7") 55764 (55cm/21.7")



**Basic Sling**  
The Basic sling is a comfortable sling without reinforcements, handles, or padding. See our entire Sling portfolio on our website.



**Handset**  
Two speed with up/down mode. Input for charging and USB interface. IPX7 water protection rating. Part no.: 71457



**Upholstered Cover**  
Protects user from hanger bar impact. Hanger bar sold separately. Part no. 50866



**Lift charger**  
Charger with holder. Easily attaches to wall. Part no. 50899-US




**Secondary Control**  
Secondary control panel on back of lift for emergency stop function. Integrated feature.



**4-point suspension**  
Enables 4-point lifting. Placed on leg hook. Part no.: 55619



**Distance Strap**  
Reinforced nylon thread and metal hooks. Extend lift strap up to 100cm/39.3". Part no. 20124-xx (specify length)

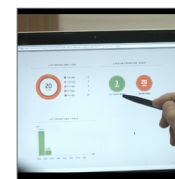

**Insights**  
Lift intelligence that facilitates cloud based storage of collected data. Integrated feature.



**Lift Trolley**  
Simplifies transport of lift between rooms. Convenient sling storage. Part no.: 70346



**Scale**  
Medical grade scale. Max. capacity 320kg/705.5lb. ISO and EN compliant. Art. no.: 50590 (RoamerS & Altair) 50591 (Heliq)

\*Spreader Bar Small: 92441 is manual.

# Technical Specifications

Description	Roamer S		
Part Number (Lift)	55300H-[Country Code]	Part Number (Hanger Bar)	55751 (45cm/17.7") 55764 (55cm/21.7")
Dimensions (W x H x L)	40cm x 16cm x 30cm 15.7" x 6.2" x 11.8"		
Safe Working Load	220kg/484lb		
Product Weight (Lift)	8.9kg / 19.6lb		
Product Weight (Hanger Bar)	1.3kg/2.8lb (55751) or 1.5kg/3.3lb (55764)		
Battery	2.3Ah, 2x12V		
Battery Charger	Input: 100-240V, AC 50-60 Hz Output: 27.8V, 0.8A DC		
IP Class	Lift: IPX4 Hand Control: IPX7		
Lifting Distance and Speed	Distance: 21.2cm/8.3" Speed: 1.8m/min - 3m/min — 70"/min - 118"/min		
Approvals	ROHS-2 approved, Tested and approved according to ISO 10535:2006 ETL Listed		
Emergency Lowering	Electronic and Mechanical		

Country Codes: AU-Australia, GB-Great Britan, EU-Europe, JP-Japan, NA-North America

