



Cableless  
Comfortable & Stylish  
High Weight Capacity!



## Our bestselling mobility aid

neXus 3, is our bestselling rollator with a wide backstrap and robust wheels, making it a perfect everyday-life-companion.



### Easy to fold

With the first-ever patented cross-folding frame, it set the standard for future walking aids design. The neXus 3 rollator features a cross-folding frame, making it exceptionally easy to fold. The soft zippered bag folds together while remaining on the rollator. It can stand independently in a folded position with the bag on the rollator. The cross-folding design also enables the rollator to be partially folded to pass through narrow spaces.



### Cableless and wide backstrap

The cable-free and anatomic brake handles increases safety and comfort. No cables catch on objects around the home. The cable-free braking system also ensures years of trouble-free use. We have had anatomy and strong

visual communication in mind when designing the brake handles. They are comfortable, and the red ball is designed for those with limited hand dexterity and allows the user to apply/disengage the brakes, using the palm of their hands. The cross-folding structure enables superior walking ergonomics with ample walking space.

The padded, contoured, and cushioned seat offers a comfortable resting spot. The backstrap can be flipped to allow for two-position seating as well as a visual aid for achieving optimal gait.


### Available in black, blue and red

The neXus 3 is available in three different sizes; neXus 3 Super Low (19"), neXus 3 Low (22"), and neXus 3 Standard (24")\*. All sizes are available in Black, Blue and Red (Burgundy) colour.

\* Measurement for sizes is floor to bottom of seat pan.



## Features

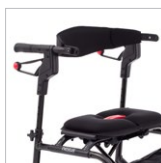


### Padded Backrest

Back comfort while seated. Can be flipped both ways and is included with the rollator.



### Foam Cushioned Seat

Offers a comfortable resting spot while walking.

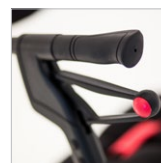

### Anatomic Brake Handles

The anatomic brake handles are easy and comfortable to use.

## Accessories




### Hard Basket

Basket ideal for shopping. Maximum load 4.5kg/10 lb. Part no. 4510.




### Tray

Practical for resting lightweight items. Maximum load 2 kg/4.4lb. Part no. 93316.

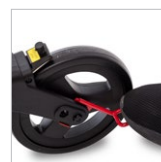


### Curb Climber

Allows the rollator to be slightly lifted onto a curb. Part no. 93315.

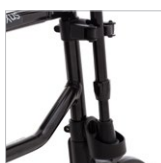


### Cane or Umbrella Holder

Simplifies transport of a cane or umbrella. Part no. 93317.

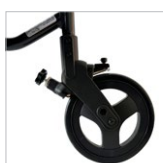

### Slow Down Brake

Enables controlled walking speed. Delivered in pairs. Part no. 93122.

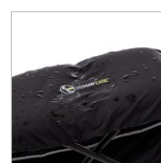

### Seat Cover

Waterproof seat cover that is easy to attach to the seat. Part no. 93313.



### Cup Holder

Attaches in seconds without tools. Fits standard size bottle and cups. Part no. 90394.



### Identification Tag Canada

Attaches to the rollator in seconds. Part no. 90392.



### Soft Zippered Bag

Soft bag that folds with the mobility aid. Maximum capacity 5kg /11lb. Part no. 93312.

## Technical specifications

Part Name	Part no.	Handle Height (Min - Max)	Seat Height	Overall Width	Width between Handles	Distance between Rear Wheels	Wheel Size	Product Weight
neXus 3 Super Low	Black 470018*	78 cm - 86 cm 30.5" - 34.1"	48 cm 19"	58 cm 23"	45 cm 17.5"	47 cm 18.5"	20 cm x 3 cm 8" x 1.2"	7.6 kg 16.9lb
	Blue 470018B							
	Red (Burgundy) 470018R							
neXus 3 Low	Black 470022	78 cm - 93 cm 30.5" - 36.5"	56 cm 22"	58 cm 23"	45 cm 17.5"	47 cm 18.5"	20 cm x 3 cm 8" x 1.2"	7.7 kg 17.1lb
	Blue 470022B							
	Red (Burgundy) 470022R							
neXus 3 Standard	Black 470024	83 cm - 101 cm 32.5" - 39.8"	61 cm 24"	58 cm 23"	45 cm 17.5"	47 cm 18.5"	20 cm x 3 cm 8" x 1.2"	7.8 kg 17.3lb
	Blue 470024B							
	Red (Burgundy) 470024R							

Maximum User Weight: 200 kg / 441 lb.

\* Article number 470018S (neXus 3 Super Low Black) can be equipped with a different handle kit to lower the handle height to 27,5"-30,5".



The neXus 3 is available in three different colours; Black, Blue, and Red [Burgundy].

**”** We improve freedom of movement and quality of care for people.